COMON



BQ80 Infant Radiant Warmer

New Comprehensive Neonatal Nursing Platform

This is by far the most powerful 4 in 1 neonatal nursing platform. The BQ80 integrates four key rescue and nursing systems to achieve one-stop operation and management. At the same time, it scientifically optimizes the workflow, helps medical staff to easily respond to urgent medical needs, effectively saves the valuable rescue time, and provides comprehensive and meticulous care for newborns.





Warm and Caring

Far infrared ceramic heating tube

Utilize Elstein far infrared ceramic heating tube, which greatly improve the radition efficiency and extend the duration.

Perfect to absorded by skin and effectively accelerate blood circulation, by which the visceral temperature and body surface temperature rise simultaneously in short time.

Double heating tube design accelerate the heating rate and ideally distribute the heat on the surface even the mattress is tilted, which provides warm and comfortable for neonatal.

Apply high-accurate double thermistor temperature sensor and increase its accuracy of temperature monitoring.

Thermal conduction gel mattress

Environmental silicon gel, with great biocompatibility Soft and comfortable, great thermo performance Anti-bacterium, easy to clean and disinfect as well as reduce bedsore





BQ80 integrates specialized neonatal monitoring system. It can monitor the vital signs and waveforms of newborns in real time during resuscitation without additional patient monitor, so that medical staff can evaluate and decide if the resuscitation measures are effective

Specialized Neonatal Monitoring System

Includes specialized ExNeo neonatal ECG technology, Adap- DSP self-adaptive blood pressure processing technology, Masimo SET neonatal blood oxygen technology, patented neonatal apnea self-rescue technology, and equipped with specialized neonatal accessories that have passed the biocompatibility test; Whether it is to ensure the accuracy of ECG, respiration, blood oxygen, blood pressure, or to treat the suffocating symptoms caused by apnea, it has obvious advantages, which provides great convenience for newborn rescue and resuscitation;





Ergonomic Design

Simplify operation, improve efficiency, and save valuable time for medical staff.

Hands-free Alarm Silence

It's able to eliminate the alarm sound just by waving your hand in front of the sensor without touching the machine; help you maintain a clean, quiet and warm nursing environment, and effectively prevent the occurrence of hospital infection.





Self-adaptive Brightness Adjustment Touch Screen

10.4-inch color LCD touch screen

Centralized display of all neonatal waveforms and data for easy viewing and operation


Ambient light sensor, automatically adjusted screen brightness

360° Visible Alarm Indicator

The alarm can be detected from every angle, whereby event can be more noticeable .



Humanized Damping Door Built-in damping components, reduce the influence of noise to neonatal

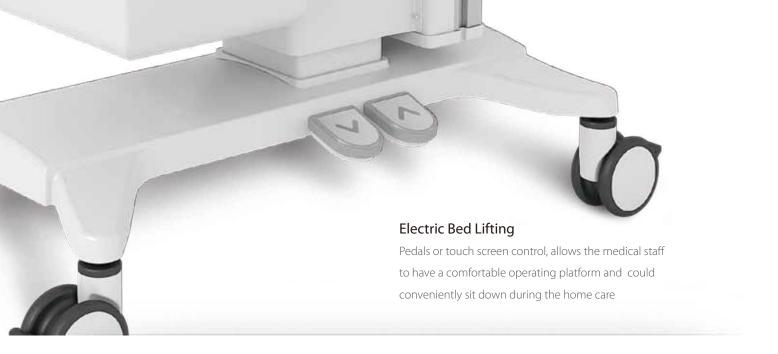



Bed Tilting

The bed can be tilted through the touch screen, support one-button angle adjustment, can be leveled with one button when weighing, accurate and fast;

Patented Tilt Angle Display

The patented tilt angle display function is accurate and visualized by built-in angle sensor,



COMEN 08

Rotary Lamp

Built-in X-ray Film Box

X-ray film box is placed under the crib, no need to move the neonatal during shooting X-rays, avoiding adverse stimulation and reducing the risk of infection.



Integrated Electronic Scale

Measuring the weight of neonatal, convenient to record and observe neonate growth trends



COMEN Share with the World

C€₁₆₃₉

Johannesburg 011 430-7000 Cape Town 021 983-1300 **East London** 043 727-1241 **Durban** 031 266-5518

Gqeberha041 363-4928 **Bloemfontein**051 448-2183

Sharecall 086 111 7736 Website www.ssemmthembu.co.za

