


Meditec M700 Series

Patient Monitors



Advanced Monitoring Parameters

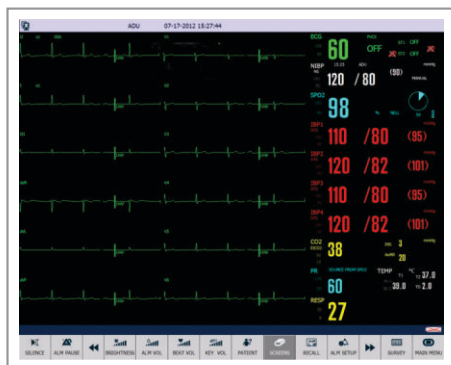

12-lead ECG waveform display

ECG-CardioTec™ 12-lead ECG Analysis

- 12-lead ECG to analyze I, II, III, aVR, aVL, aVF, V1-V6 and display these waveforms on the same screen
- Two separate ECG cable to get accurate 12-lead ECG analysis



SpO₂ Masimo


- The first signal processing and sensor technology capable of accurately monitoring physiology during motion artifact and low perfusion

Gold Standard of SpO₂

- Identifies the venous blood and other noise, isolates them, and using proprietary techniques, extracts the arterial signal



Mainstream CO₂ & Sidestream CO₂

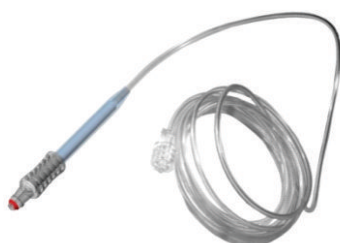

Nomoline™




Adult/Pediatric
and Infant airway adapters

Advanced Plug-In Modules

Multi-GAS - Masimo Sidestream & Mainstream Para-magnetic O₂ sensors

- Long life cycle
- Saving ongoing hospital cost

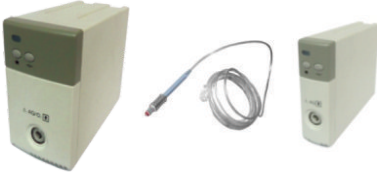


Support 4-Channel IBP measurement

- Abbott / Edwards/ BD / Utah IBP accessories can be selected to monitoring ABP, PAP, CVP, ICP etc



Electroencephalography EEG measurement

- Provides realtime waveform from 4 channels / 8 channels / 16 channels
- EEG module is designed for continuous, real-time monitoring of adult, paediatric and neonatal patients in anaesthesia and intermediate/critical care environments.

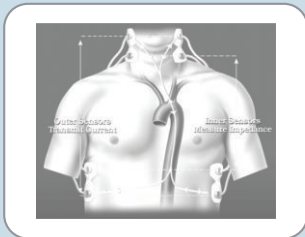


ICG - US Bloz® Non-invasive Hemodynamic

- Impedance cardiograph for non-invasive continuous hemodynamic monitoring
- Micro-signal transmission through disposable electrode
- Measurement and display of ICG
- Blood volume and Blood Flow Velocity varies with heartbeat. DISQ® technology processes impedance signal variation
- Variation of impedance applies to non-invasive Z MARC algorithm for acquiring SV, CO, SVR, Contractility and TFC etc.



M797 & M757 modular monitors with option of transport monitor and separate modules for all vital signs


Calculation Screens for Clinical use

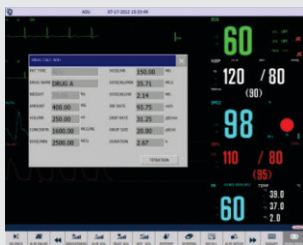

Drug




Oxygenation

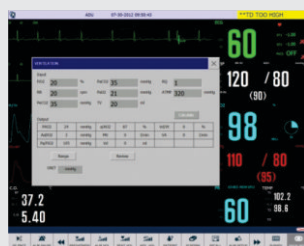

Ventilation



Kidney

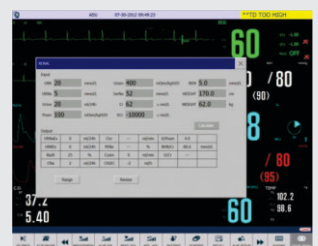

Hemodynamic



Module Expansion Rack
(M797, M757)





	M797	M777	M771	
TFT Screen Size	17"	15"	19"	
Touch Screen Operation	✓	✓	✓	
Type	Modular	Upgradable	Modular	
ECG Leads	3/5	3/5	3/5	
12 Lead ECG Monitoring & Analysis	(Optional)	(Optional)	(Optional)	
ECG Review - Arrhythmia & ST Segment Analysis	✓	✓	✓	
SpO ₂	Self/ Masimo	Self/ Masimo	Self/ Masimo	
Temp Channels	2	2	2	
IBP	Up to 4 (optional)	Up to 4 (optional)	Up to 4 (optional)	
EtCO ₂	(Optional)	(Optional)	(Optional)	
Anaesthetic Gas Module (CO ₂ , N ₂ O, HAL, ISO, ENF, DES, SEV, MAC)	(Optional)	(Optional)	(Optional)	
Paramagnetic Oxygen Monitoring	(Optional)	(Optional)	(Optional)	
Cardiac Output	(Optional)	(Optional)	(Optional)	
Depth of Anaesthesia / BIS	(Optional)	(Optional)	(Optional)	
Respiratory Mechanics (RM)	(Optional)	x	(Optional)	
SV02	(Optional)	x	(Optional)	
ICG - Impedance Cardiography	(Optional)	x	(Optional)	
Neuromuscular Transmission (NMT)	(Optional)	x	(Optional)	
EEG	(Optional)	(Optional)	(Optional)	
Oxy CRG	✓	✓	✓	
Waveforms	Up to 12	Up to 12	Up to 12	
Waveform Recall Function	✓	✓	✓	
Drug Calculator	✓	✓	✓	
Graphics & Tabular Trends	120 hrs	120 hrs	120 hrs	
Defibrillation Interface & Cautery Protection	✓	✓	✓	
Pacemaker Detection	✓	✓	✓	
Built-in Li-ion Battery Backup	2 hrs	4 hrs	2 hrs	
Data Storage and Review	4 USB & SD Card	2 USB & SD Card	4 USB & SD Card	
HL7 Compliant	✓	✓	✓	
Nurse Call Compatible	✓	✓	✓	
Connectivity with Central Monitoring System	Wired/Wireless connectivity with Meditec CMS	Wired/Wireless connectivity with Meditec CMS	Wired/Wireless connectivity with Meditec CMS	
In-built Thermal Printer	(Optional)	(Optional)	(Optional)	
Alarm System (Audible & Visual)	3 Level audio/ Visual alarm	3 Level audio/ Visual alarm	3 Level audio/ Visual alarm	
Net Weight	11 Kg	4.8 Kg	6 Kg	
Module Expansion Rack	(Optional)	x	(Optional)	
Wall Mount	(Optional)	(Optional)	Standard	
Mobile Trolley	(Optional)	(Optional)	(Optional)	
Bed Rail Mount	(Optional)	(Optional)	(Optional)	
Option for Transport Monitor	✓	x	x	
Bar Code Scanner	(Optional)	(Optional)	(Optional)	



M757	M752	M751	M747	M707
12.1"	12.1"	12.1"	12.1"	4.3"
✓	✓	✓	(Optional)	✓
Modular	Modular	Upgradable	Upgradable	Pre-Configured
3/5	3/5	3/5	3/5	3/5
(Optional)	(Optional)	(Optional)	(Optional)	x
✓	✓	✓	✓	✓
Self/ Masimo	Self/ Masimo	Self/ Masimo	Self/ Masimo	Self/ Masimo
2	2	2	2	2
Up to 4 (optional)	Up to 2 (optional)	Up to 4 (optional)	Up to 2 (optional)	x
(Optional)	(Optional)	(Optional)	(Optional)	(Optional)
(Optional)	Any One	(Optional)	x	x
(Optional)		(Optional)	x	x
(Optional)		(Optional)	x	x
(Optional)		(Optional)	x	x
(Optional)	x	x	x	x
(Optional)	x	x	x	x
(Optional)	x	x	x	x
(Optional)	x	x	x	x
(Optional)	x	(Optional)	x	x
✓	✓	✓	x	x
Up to 12	Up to 8	Up to 12	Up to 12	Up to 2
✓	✓	✓	✓	x
✓	✓	✓	✓	x
120 hrs	120 hrs	120 hrs	120 hrs	48 hrs
✓	✓	✓	✓	✓
✓	✓	✓	✓	✓
4 hrs	4 hrs	4 hrs	5 hrs	8 hrs
4 USB & SD Card	2 USB	2 USB & SD Card	1 USB	x
✓	✓	✓	✓	x
✓	x	✓	✓	x
Wired/Wireless connectivity with Meditec CMS	Wired/Wireless connectivity with Meditec CMS	Wired/Wireless connectivity with Meditec CMS	Wired/Wireless connectivity with Meditec CMS	Wired/Wireless connectivity with Meditec CMS
(Optional)	(Optional)	(Optional)	(Optional)	x
3 Level audio/ Visual alarm	3 Level audio/ Visual alarm	3 Level audio/ Visual alarm	3 Level audio/ Visual alarm	3 Level audio/ Visual alarm
5.9 Kg	3.4 Kg	4.8 Kg	3 Kg	Light Weight 1.2 Kg
(Optional)	(Optional)	x	x	x
(Optional)	(Optional)	(Optional)	(Optional)	x
(Optional)	(Optional)	(Optional)	(Optional)	x
(Optional)	(Optional)	(Optional)	(Optional)	(Optional)
✓	x	x	x	x
(Optional)	(Optional)	(Optional)	(Optional)	(Optional)

Technical Specifications - Standard Parameters

Display	
General Parameters	Suitable Bedside monitors for Adult/Paediatric & neonatal patients User selectable priority alarms - Red, Yellow, Orange color coded
Sweep Speed	6.25mm/s, 12.5mm/s, 25mm/s, 50mm/s
Various working interface selectable	Standard Monitoring Display, Large Font Intensive Care Display, Trend Graph/ Monitoring Co-Display, Bed to Bed view Display (optional), OxyCRG Dynamic View Display, Drug Dose Calculation Interface
Trends	Tabular & Graphical trends
Events	ECG, SpO ₂ , NIBP with waveform & review

ECG	
Lead Mode	3 Leads (I, II, III), 5 Leads (I, II, III, aVR, aVL, aVF, V) 12 Leads ((I, II, III, aVR, aVL, aVF, V1-V6)
Waveform	Up to 12 waveforms display on screen
4 selectable plus	0.25, 0.5, 1, 2 cm/mV
PACE Detection	Automatic
Defibrillation energy distributary	
Input voltage range	±8mV, polarized voltage: ±500mV
CMRR Diagnosis	> 100 dB (no 50Hz/60Hz software wave trap)
Monitor	> 110 dB (has 50Hz/60Hz software wave trap)
Surgery	> 100dB (50Hz/60Hz software wave trap)
System noise	<30μV p-p (RTI)

HR Measuring and Alarm Range	
Adult/Paed	15 bpm – 300bpm
Neo	15 bpm – 350bpm
Accuracy	±1% or ±1bpm, which great
Resolution	1 bpm
Sensitivity	> 200 μV P-P
Differential Input Impedance	≥ 5 MΩ
Electrode offset potential	300mVd.c. – 600mVd.c.
Patient leakage current	< 10 uA
ECG Signal Range	6 mV (Vp-p).
ST Segment Monitoring Range	
Measure and Alarm	-2.0 – +2.0 mV

Respiration			
Method	Impedance between R-F (RA-LL), R-L (RA-LA)		
Respiration Rate Measuring and Alarm Range:			
Adult	0 rpm – 120 rpm	Neo/Ped	0 rpm – 180 rpm
Resolution	1 rpm	Accuracy	2 rpm
Gain Selection	0.25,0.5, 1, 2, 3, 4, 5		

Technical Specifications - Optional Parameters

IBP			
Channel	Up to 4 channel		
Label	ART, PA, CVP, RAP, LAP, ICP, P1, P2		
Static Pressure Measuring Range	(-)50 – +300 mmHg (up to 350 mmHg)		
Static Pressure Accuracy	±2% or ±1mmHg whichever is higher		
Dynamic Pressure Measuring Range	(-)50 – +300 mmHg		
Dynamic Pressure Accuracy	2 % or 1mmHg which bigger		
Measuring and Alarm Range			
ART	0 – 300 (mmHg)	PA	6 – 120 (mmHg)
CVP/RAP/LAP/ICP	(-)10 – 40 (mmHg)	P1/P2	(-)50 – 300 (mmHg)
Resolution	1 (mmHg)	Zero range	200 mmHg

NIBP	
Method	Oscillometric
Mode	Manual, Auto, Continuous
Measuring Interval in AUTO mode	1/2/3/4/5/10/15/30/60/90/120/240/480 Min
Continuous	5min, interval 5s
Measuring Type	Systolic Pressure, Diastolic Pressure, Mean Pressure
Cuff Pressure measuring Range	0–300 mmHg
Pressure Resolution	1 mmHg
Mean error	3 mmHg
Maximum Standard deviation	3 mmHg
Over pressure protection	Dual Over pressure protection
Measured unit	mmHg/kPa selectable
Measurement types	Systolic, Diastolic, Mean Systolic : 40-300 mmHg Diastolic : 10-215 mmHg Mean : 20-235 mmHg
Pressure range for Paediatrics	Systolic : 20-235 mmHg Diastolic : 10-150 mmHg Mean : 20-165 mmHg
Pressure range for Neonates	Systolic : 10-235 mmHg Diastolic : 10-100 mmHg Mean : 20-110 mmHg

SpO ₂	
Method	Self / Masimo
Measuring Range	1 – 100%
Alarm Range	1 – 100%
Resolution	1 %
Accuracy with Adult / Paediatric with Low-perfusion	1 digit (70%–100% SpO ₂) Undefined (0–70% SpO ₂)
Neonate	1 digit (70%–100% SpO ₂) Undefined (0–70% SpO ₂)
Pulse Rate Measuring and Alarm Range	20 – 300bpm

Temp			
Channel	2	Measuring Range	0 – 50°C
Sensor type	YSI (B Series) and CF-FI	Resolution	0.1°C
Accuracy	0.1°C (25°C–45°C), 0.2°C (0–25°C, 45° – 50°C)		

CO ₂			
Channel		Infra-red Absorption Technique	
Measuring Mode		Sidestream or Mainstream	
Measuring Range			
CO ₂	0 – 150 mmHg	AwRR	2 – 150 rpm
Unit	mmHg, %, Kpa	Resolution	1 mmHg, 0.1%, 0.1Kpa
CO ₂	1mmHg	InCO ₂	1mmHg
Accuracy - AwRR		±1 rpm	
Accuracy - CO ₂		2 mmHg, 0 – 40 mmHg Reading 5%, 41 – 70 mmHg	
Alarm Delayed - AwRR		10 – 40 seconds	
Response time		<3 seconds, includes transport time, risetime	
Resolution		1%	

Technical Specifications - Optional Parameters

Cardiac Output (CO)

Method	Thermodilution Technique		
Measuring range - CO	0.1 – 20 L/min		
TB	23 – 43°C	TI	(-)1 – 27°C
Resolution - CO	0.1 L/min	TB, TI	0.1°C
Accuracy - CO	5% or 0.2 L/min	TB	0.1°C
TI	0.1°C		
Output parameters	CO, Hemodynamic Calculation		
Alarm range	23 – 43°C		

Neuromuscular Transmission (NMT)

Microprocessor-controlled	
Stimulator Mode	TOF, TOFS, PTC, 1Hz Twitch, 0.1Hz Twitch, DBS 3.3 and 3.2 (double burst), Tetanic Stimulation (burst), 5s-50Hz or 100Hz Output (accuracy ±5% of full scale value)
Surface Electrodes	Constant current, 0-60mA(0-12/18μC) up to 5Kohm, Monophasic, 200μs or 300μs pulse width
Needle Electrodes	Constant current, 0-6mA (0-0.24μC) up to 5KOhm. Monophasic, 40μs pulse width
Acceleration Transducer	Accuracy ±5% of full scale value
Temp. Sensor	Range 20.0-41.5°C (accuracy ±5°C)

Multigas (AG)

Technology	Side Stream or Main Stream		
Gas	CO ₂ , N ₂ O, Des, Iso, Hal, Sevo, Enf, O ₂		
Measuring Range			
CO ₂	0 – 10%	O ₂	0 – 100%
N ₂ O	0 – 100%	AwRR	2 – 100 rpm
Halothane	0 – 5%	Isoflurane	0 – 5%
Enflurane	0 – 5%	Sevoflurane	0 – 8%
Desflurane	0 – 18%		
Other	Up to 3 waveforms display, Agent mixture detection MAC value display		

Depth of Anaesthesia/ BIS

Bispectral Index (BIS)	0 – 100
Electromyographic Strength (EMG)	30 – 80 dB (where 1μV ² = 40dB)
Signal Quality Index (SQI)	0 – 100%
Suppression Ratio (SR)	0 – 100%
Spectral Edge Frequency (SEF)	0.5 – 30.0Hz
Total Power	40 – 100dB (where 1μV ² = 40dB)
Burst Count	0 – 30 (with an Extend Sensor only)
EEG (Continuous, single channel real-time waveform(Standard) / 4 channel (Optional)	
Noise (EEG Waveform)	< 0.3 μV RMS (2.0 μV peak-to-peak) : 0.25 Hz – 50 Hz
EEG Bandwidth	0.25 Hz – 100 Hz (-3 dB) ± 10%

Recorder

Built-in, direct thermal pixel array recorder up to 3 channel printing	
1, 2, 3 channel selectable	
Print speed	25mm/s, 50mm/s
Paper width	50mm

Power Supply

Input Voltage	AC 100V-240V
Input Frequency	50/60Hz

Environmental Conditions

Storage	Temp	(-) 10°C – 40°C
	Humidity	15%–80%
Operational	Temp	0°C – 40°C
	Humidity	15%–80%
Atmospheric Pressure	700kPa–1060kPa	

Technical Specifications - Central Monitoring System

Standard Configuration

Processor	Dual Core
Display	17"/19"/21"/24" TFT screen
Print	HP printer

Features

Resolution	1280 × 1024
Language	English
Waveform	ECG, RESP, SpO ₂ , NIBP, EtCO ₂ , etc.
Parameter	HR, RR, ST, SpO ₂ , PR, NIBP, TEMP, etc.

Alarm

Alarms	High and low limits alarm
	Arrhythmia alarm
	Audible and visual alarm

Bedside Monitors

Support maximum 128 bedside monitors
Display maximum 32 patient's information
Key monitoring interface, display 6 waveforms
Intellectualized Bi-directional Communication

Data

720-hour full disclosure waveform review
1240-hour trend graphic review
10,000 alarm record review
50,000 historical patient records

Recorder

External thermal array
Real-time recording, review recording and historical
Reports recording

Connection Mode


Complete solution for LAN, wireless LAN, and mixture networking

Power Supply

Input Voltage	AC 100V-240V
Input Frequency	50Hz ± 1Hz

Central Monitoring System

Bi-directional communication with Meditec Central Monitoring Station by wired or wireless connection

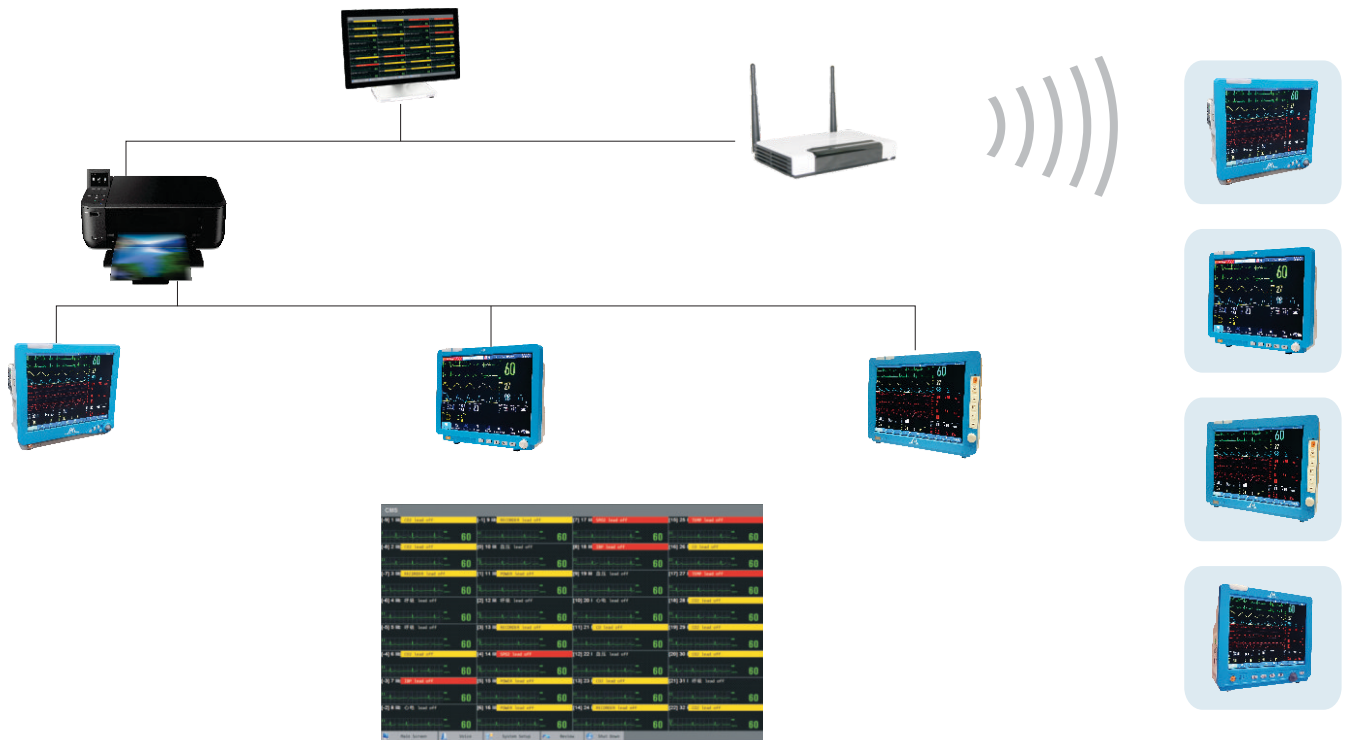

32 Beds display



Wave review and print preview



Large font display

- **Standard Network Interface:** TCP/IP central monitoring system, be able to observe full vital signs information of the patient
- **Full Disclosure Waveform Observation:** Large screen display, be able to observe 16 bedside monitors waveforms simultaneously. Selectable Dual-screen display.
- **Bi-directional Communication:** Unique Bi-directional communication function realize the mutual control of central monitor and bedside monitors
- **Large Storage Capacity:** 720-hour full disclosure waveform review, 1240-hour trend graph review, 1000 alarm records review. 50,000 patient historical data. All the information can be stored in disk.
- **Data File Review:** Large storage capacity of full monitoring information for easy review, analysis and statistics
- **Dual Screen:** Support various data display methods and various windows layout. It can display the waveform and numeric information of ECG, SpO2, RESP, NIBP, IBP, CO2, etc. The key bedside monitoring function provides more detail information to the doctor.
- **Operating System:** User-friendly central monitoring software built on Windows platform. Advanced information integration technology to realize the concentration management of all parameters. Full digital transmission to ensure the communication accurate and reliable.
- **Stable Internet Connection:** Standard configuration for 8 bedside monitors, support maximum 128 units



Meditec International England Limited

8, Pinner View, Harrow, Middx,
HA1 4QA, United Kingdom
Tel: +44 203 6423808

E-mail : sales@meditecengland.co.uk
Website : www.meditecengland.co.uk

"In view of continued product improvement Meditec International England Limited reserves the right to make changes in the equipment which may affect the information contained in this leaflet without giving any prior notice. Pictures are representative, accessories are optional and may vary at the time of supply."